


Visualize the whole heart with DualFlow,  
our latest innovation for the Stellant D

*"If my practice were beginning coronary CTA today, I would definitely invest in an injector that can inject both contrast and saline."*



## Stellant® CT Injection Systems


**No Saline Flush**

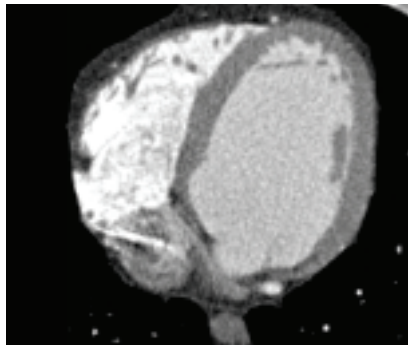

**with Saline Flush**



**with Simultaneous Injection**

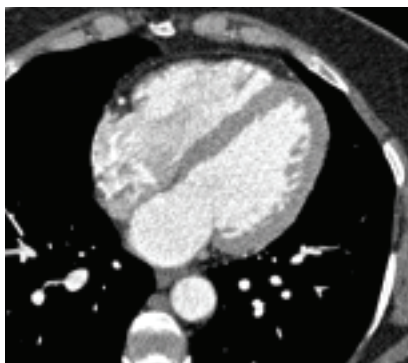


### **Variable Motion Ratios**

Stellant D with DualFlow option is more than a saline flush after a contrast bolus. Now you can inject both contrast and saline at the same time. The key is the simultaneous injection capability of the DualFlow option. DualFlow enables variable ratios of plunger motion from the contrast and saline syringes simultaneously. With the proper ratio, left and right heart ventricles can be illuminated uniformly for improved image quality. <sup>1</sup>

### **With DualFlow you can:**

- Match attenuation levels in both ventricles
- Minimize artifacts
- Visualize right coronary arteries and right ventricle in a single study

MEDRAD CT...Enhancing Your CT Practice.

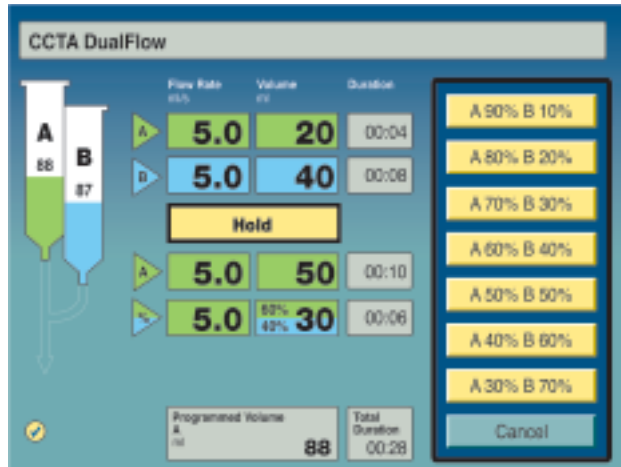

## Simple Upgrade to DualFlow

It's easy to unlock Stellant D with DualFlow option. Each Stellant D already has the capability for DualFlow simultaneous injection.\* To upgrade, just contact your MEDRAD Customer Support Representative for the code to unlock the DualFlow option on your Stellant D. Enter the code and begin using the DualFlow option for your cardiac CT studies.

\*Some systems may require an upgrade by MEDRAD Service Engineers

<sup>1</sup>Sablaylorles, J-L, "Cardiac CT: Experience from Daily Practice", Advanced CT, Aug 2004: 4-10.



Catalog Number: SCT 221  
Stellant D Injection System with Pedestal Mount and Dual Flow

Catalog Number: SCT 222  
Stellant D Injection System with Pedestal Mount and Dual Flow

Catalog Number: UFK D2DF  
DualFlow Upgrade Kit

MEDRAD reserves the right to modify the specifications and features described herein, or discontinue manufacture of the product described at any time without prior notice or obligation. Please contact your authorized MEDRAD representative for the most current information.

MEDRAD, and Stellant are federally registered trademarks, and Performance. For life. is a trademark of MEDRAD, INC. U.S.A.



© 2005 MEDRAD, INC. All Rights Reserved.

STL SI NA 203216 Rev B



the Baldrige Award recognizes exemplary achievements in seven areas: leadership, strategic planning, customer and market focus, measurement, analysis and knowledge management, human resource focus, process management and business results. The U.S. Commerce Department's National Institute of Standards and Technology (NIST) manages the Baldrige National Quality Program in close cooperation with the private sector. The award is traditionally presented by the President of the United States in a special Washington, D.C. ceremony.

The Malcolm Baldrige National Quality Award is our nation's premier award for performance excellence and quality achievement. Established in 1987 and named after the former Secretary of Commerce,

**MEDRAD, INC.**  
One Medrad Drive  
Indianola, PA 15051-0780 USA  
412-767-2400

**Customer Service/Orders**  
1-800-MEDRAD-1  
(1-800-633-7231)

**Customer Service FAX**  
(412) 767-4120

**International FAX**  
(412) 767-4128

**For Service Repairs**  
1-800-MEDRAD-7  
(1-800-633-7237)

**MEDRAD Service FAX**  
(412) 767-4126

**www.medrad.com/stellant**  
Email address for more information:  
info@medrad.com

